


INFORMATIVE NOTE
Name: REFLEXX 72
Date: 31/05/2023 Rev: 09
MANUFACTURER

 REFLEXX SpA
 Via Passeri, 2 46019 VIADANA (MN)
 Tel.: + 39 0375 758891 Fax + 39 0375 464504 info@reflexx.com
www.reflexx.com

Imballo - Pack

100		10	
Taglia Size	Articolo Article	Taglia Size	Cartone Carton
R72 S	8032891632598	S	8032891637593
R72 M	8032891632604	M	8032891637659
R72 L	8032891632611	L	8032891637616
R72 XL	8032891632628	XL	8032891637623

DESCRIPTION: ■ disposable glove powder free, not sterile ■ **Material:** nitrile ■ **Shape:** ambidextrous; five fingers, beaded cuff ■ **Colour:** light blue ■ **Feature:** internally chlorinated; textured surface and textured fingertips ■ **Packing:** 100 pieces/box; 10 boxes/carton ■ **ISO 2859:** AQL 1.5 (level G1) freedom from holes; AQL 2.5 (level G1) major defects; AQL 4 (level G1) minor defects

C CONFORMITY AND PERFORMANCES
Personal Protective Equipment cat III


0465 ANCI_CIMAC, Via Aguzzafame, 60/b 27029 Vigevano (PV)

Reg EU 2016/425 - EN ISO 374-1:2016 TYPE C - EN ISO 374-5:2016 - EN 374-2:2014 - EN 420:2003 + A1 2009

Medical Device

Reg EU 2017/745

(CND T01020204 NON SURGICAL NITRILE GLOVES)

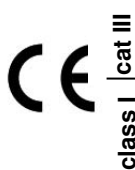


Reg CE 1935/2004 – Reg EU 10/2011 - DM 21/03/73. To such use carefully read the declaration of suitability for food contact (DICA).

DIMENSIONS		
Cod./Size	Length min(mm)	Width (±5 mm)
R72/ S (7)	240	85
R72/ M (8)	240	95
R72/ L (9)	240	105
R72/ XL (10)	240	115

Thickness mm	Palm (+/- 0.02)	0.07
--------------	-----------------	------

EN 455-2:2009/A1	Tensile strength	≥6 N
------------------	------------------	------

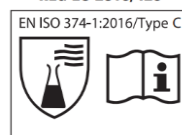REG EU 2016/425

EN ISO 374-5:2016


EN ISO 374-1:2016				
Chemical compound :	EN 16523-1:2015 TYPE C		EN 374-4:2013	
	COD	Performance level	Breakthrough time	% degradation
Sodium hydroxide 40%	K	3	>60'	7,1%

RISKS FROM WHICH THE DPI IS DESIGNED TO PROTECT:

Chemical compound :	Protection
K Sodium hydroxide 40 %	CAS:1310-73-2 Protection from use of inorganic bases

RECOMENDATIONS OF USE Multi-purpose glove to protect the hand and/or the handled product, can also be used for general uses ■ The glove is suitable for contact with pure chemical products for a period of less than the permeation time indicated in the chemical resistance table ■ Results obtained under laboratory conditions may not reflect actual working time **WARNINGS** Before using the product, try it in the real condition of use ■ The glove is a device subject to wear, so frequent replacement is recommended and in all cases where imperfections are observed ■ Wear gloves with dry, clean hands ■ In rare cases where transient hypersensitive skin reactions may occur, discontinue use and consult a doctor ■ The user must inform the manufacturer and the relative authority of any incidents occurring during the use of the protective equipment **STORAGE AND DISPOSAL** Store the gloves in their original packaging in a dry place. ■ Do not expose to direct sunlight and heat sources. ■ The expiry date is 3 years from the date of



I N S A F E H A N D S

manufacture indicated on the sales package, and is valid for the product properly stored ■ Gloves must be disposed of according to local regulations ■ Product Declaration of Conformity on www.reflexx.com

REFLEX S.P.A. UNIPERSONALE Via Passeri, 2 - 46019 Viadana (MN) | TEL. + 39 0375 758891

Cap.Soc. 1.200.000 € (i.v.) | PIVA 02085450209 | M5UXCR1 | R.E.A. 223166 | info@reflexx.com www.reflexx.com



5010014617